

Product Data Sheet

Anti-NF-kappaB p105/p50 Antibody

Catalog # Source Reactivity Applications

CPA4479 Rabbit H, M, R, B, C, D, P WB, IH, IF/IC

Description Rabbit polyclonal antibody to NF-kappaB p105/p50

Immunogen KLH-conjugated synthetic peptide encompassing a sequence within the center

region of human NF-kappaB p105/p50. The exact sequence is proprietary.

Purification The antibody was purified by immunogen affinity chromatography.

Specificity Recognizes endogenous levels of NF-kappaB p105/p50 protein.

Clonality Polyclonal

Conjugation

Form Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol,

and 0.01% sodium azide.

Dilution WB (1/500 - 1/1000), IH (1/50 - 1/100), IF/IC (1/50 - 1/200)

Gene Symbol NFKB1

Alternative Names Nuclear factor NF-kappa-B p105 subunit; DNA-binding factor KBF1; EBP-1; Nuclear

factor of kappa light polypeptide gene enhancer in B-cells 1

Entrez Gene 4790 (Human); 18033 (Mouse)

SwissProt P19838 (Human); P25799 (Mouse); Q63369 (Rat)

Storage/Stability Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid

freeze/thaw cycles.


Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC- Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish

COHESION BIOSCIENCES LIMITED


WEB ORDER SUPPORT CUSTOM
www.cohesionbio.com order@cohesionbio.com techsupport@cohesionbio.com custom@cohesionbio.com




Product Data Sheet

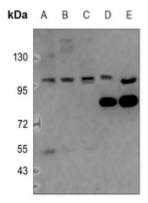

Western blot analysis of NF-kappaB p105/p50 expression in MCF7 (A), U2OS (B), A375 (C), mouse muscle (D), mouse liver (E) whole cell lysates. (Predicted band size: 85; 105 kD; Observed band size: 105; 50 kD)

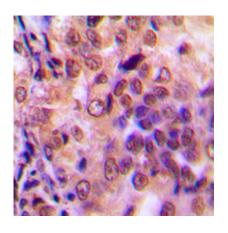

Immunohistochemical analysis of NF-kappaB p105/p50 staining in human prostate cancer formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

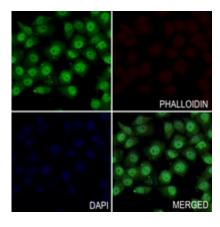

Immunofluorescent analysis of NF-kappaB p105/p50 staining in LO2 cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a hidified chamber. Cells were washed with PBST and incubated with a AF488-conjugated secondary antibody (green) in PBS at room temperature in the dark. Phalloidin - AF594 was used to stain Actin filaments (red). DAPI was used to stain the cell nuclei (blue).

Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC- Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish

COHESION BIOSCIENCES LIMITED