

### **Product Data Sheet**

# **Anti-Peroxiredoxin 5 Antibody**

Catalog # Source Reactivity Applications

CQA3174 Rabbit H, M, R WB, IF/IC

**Description** Rabbit polyclonal antibody to Peroxiredoxin 5

Immunogen Recombinant full length protein of human Peroxiredoxin 5

**Purification** The antibody was purified by immunogen affinity chromatography.

**Specificity** Recognizes endogenous levels of Peroxiredoxin 5 protein.

**Clonality** Polyclonal

Conjugation

Form Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol,

and 0.01% sodium azide.

**Dilution** WB (1/500 - 1/1000), IF/IC (1/50 - 1/200)

Gene Symbol PRDX5

Alternative Names ACR1; Peroxiredoxin-5 mitochondrial; Alu corepressor 1; Antioxidant enzyme B166;

AOEB166; Liver tissue 2D-page spot 71B; PLP; Peroxiredoxin V; Prx-V; Peroxisomal

antioxidant enzyme; TPx type VI; Thioredoxin peroxidase PMP20; Thioredoxin

reductase

Entrez Gene 25824 (Human); 54683 (Mouse); 113898 (Rat)

SwissProt P30044 (Human); P99029 (Mouse); Q9R063 (Rat)

Storage/Stability Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid

freeze/thaw cycles.


Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC- Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish

#### **COHESION BIOSCIENCES LIMITED**


WEB ORDER SUPPORT CUSTOM
www.cohesionbio.com order@cohesionbio.com techsupport@cohesionbio.com custom@cohesionbio.com



## **Product Data Sheet**

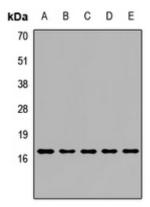

Western blot analysis of Peroxiredoxin 5 expression in Hela (A), HepG2 (B), mouse liver (C), mouse heart (D), rat brain (E) whole cell lysates. (Predicted band size: 12; 17; 22 kD; Observed band size: 17 kD)

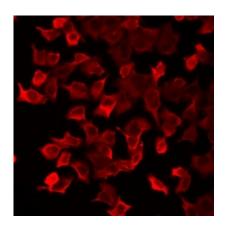

Immunofluorescent analysis of Peroxiredoxin 5 staining in U2OS cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a humidified chamber. Cells were washed with PBST and incubated with a DyLight 594-conjugated secondary antibody (red) in PBS at room temperature in the dark.

Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC- Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish

#### **COHESION BIOSCIENCES LIMITED**